

Selecting your Pacta[®]-LC

1. Consider Corporate Health & Safety


Lifts of 5 liter bottles above shoulder height should be avoided.

Top edge of the LC-system is ≥ 130 cm?

AND solvent bottles are ≥ 5 L?

→ Choose Pacta FLY.

To avoid bent back during service, bottom of LC-system should be at a height of 75-90 cm. → Choose height adjustable Pacta.

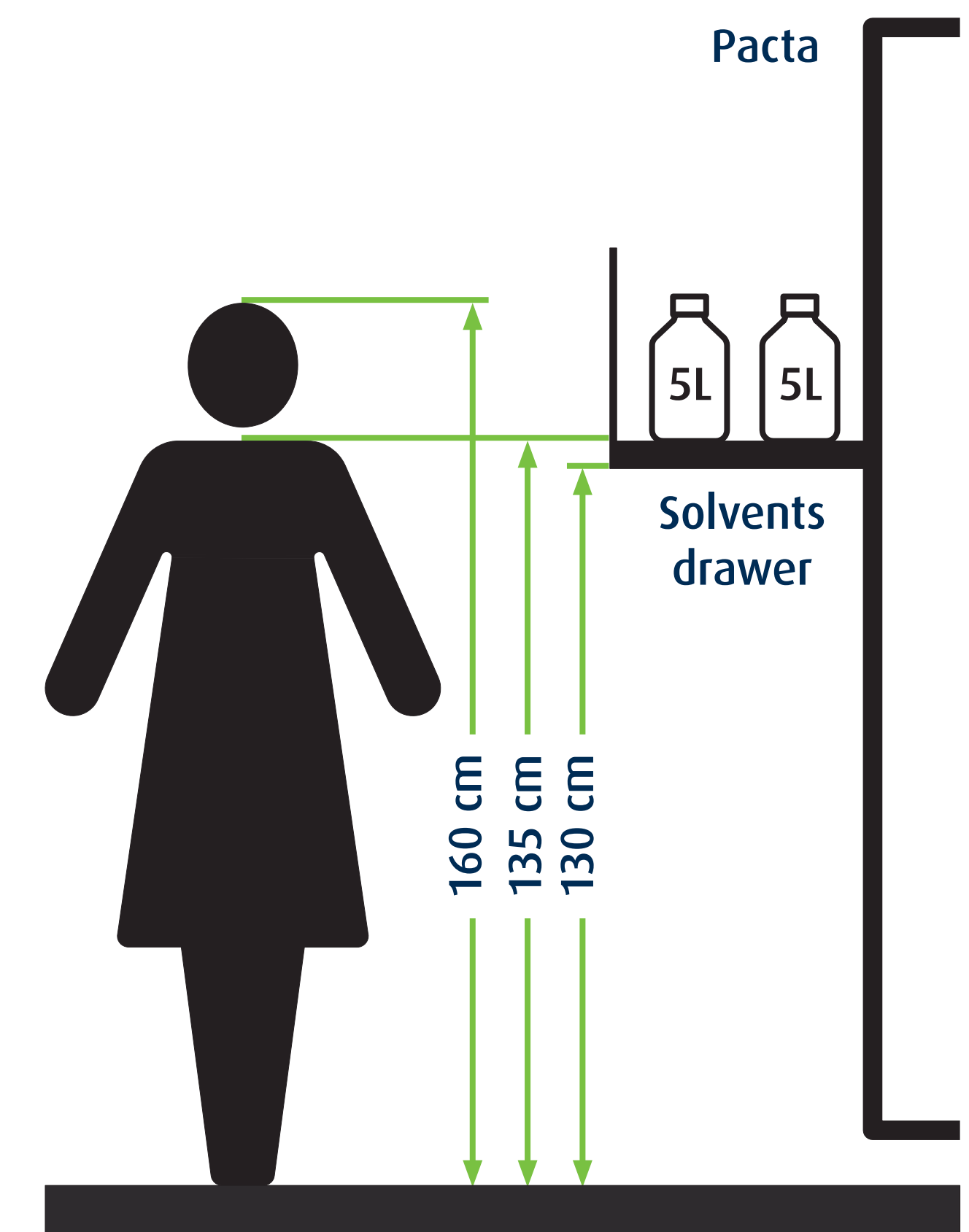


2. Determine Pacta[®] width

Pacta width is determined by LC-system width and Preparation Guide respect distance:

$$W_{\text{Pacta}} \geq W_{\text{LC}} + W_{\text{PG respect}}$$

Pacta width \geq LC system width + Preparation Guide respect distance

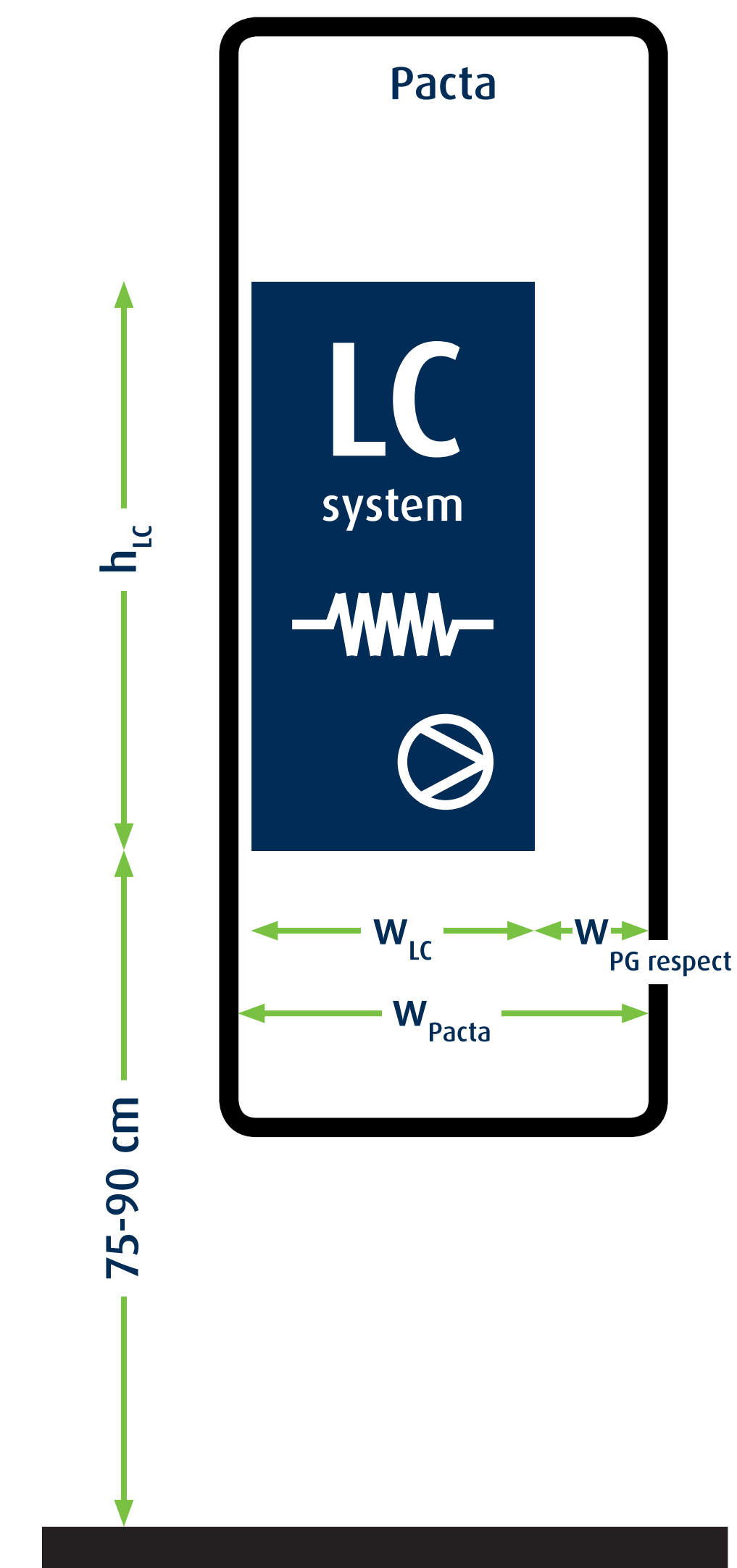

3. Determine Pacta[®] options

LC service access and LC options determine Pacta options:

Service access at the top of the LC-system?
→ Choose the Pacta detector shelf option.

Oven and service access from the side?
→ Choose the Pacta oven-swing-out option.

